



PostureFit Lift Recline Chair

Configurable Pressure Care Seating Solution for unique and changing needs



Australia
1300 133 120 | aspirecare.com.au

United Kingdom
01622 541235 | aidacare.co.uk

Personalise your comfort, optimise your posture

PostureFit Configurable Lift Recline Chairs offer adjustable and adaptable seating options designed to provide unparalleled postural support and comfort.

- 

Upright Position
- 

Lift to Stand
- 

Power Headrest
- 

Back Recline
- 

Leg Raise
- 

Full Recline



Targeted Postural Support and Interchangeable Backrests

All backrest cushions have adjustable filling and support wedges are also available to meet the unique postural requirements of the user.

Standard Backrests



Waterfall



Easily adjustable to
Kyphosis configuration

Optional Backrests



Classic



Advanced Immersion

Vertical Lifting Mechanism

First 60mm of range is vertical before commencing the tilt-to-stand motion. This creates variable seat height for taller users.

Infinite Positioning

Adjustable Headrest

Electrically adjustable via the handset, the headrest cushion can be tilted forward and backward for optimal support. Perfect neck support with the touch of a button on the handset.

Zero Gravity Positioning

Creates a feeling of weightlessness and redistributes pressure off the lower back for pain management.

Oedema Management

Leg elevation above the heart promotes drainage of lower limb swelling, aids circulation, and reduces fatigue.

No more kneeling down or manual adjustments.

Integrated Footrest Extension

Suits most user heights without the need to alter the legrest or height of the chair. The leg extension can also be easily removed for users with specific heel care needs.

Pressure Relief Seating



Hybrid Cushion Technology: Unique foam-in-cell air cushion automatically redistributes pressure and maximises user comfort.

Tool-free Seating Adjustment

Therapists or carers can quickly and easily make adjustments without the need for any tools.

Seat Depth

Backrest can be quickly and easily adjusted to achieve two seat depths with a 50mm range.

Seat Width

Armrests are interchangeable for left/right side, simply swap sides to achieve two seat widths with a 50mm range.

Move with Ease

Effortless Mobility and Safety

Large directional castors are engaged by activating a rear carer foot pedal. These castors are not in contact with the floor when in a seated position, ensuring stability and preventing accidental tipping or sliding.

Transport mode

For safety, all electrical functions are inactivated while castors are engaged to prevent the chair being activated while mobile.

Battery Backup

Ensures chair can be used in the event of a power outage or if moved between rooms during the day.




After-Market Cushion Integration:


Pre-cut seat structure accommodates third-party cushions such as ROHO without compromising performance or infection control.

Medical Grade Foam: High density foam allows for prolonged periods of sitting.

Pressure Mapping

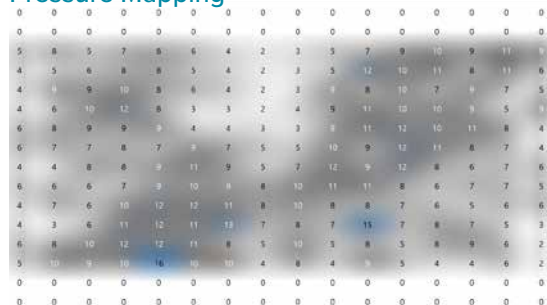

Patient A 64kg | 165cm

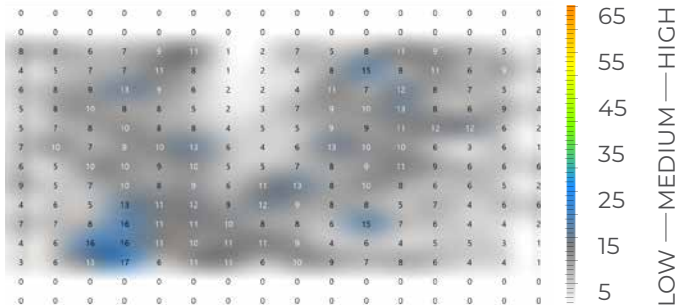

Patient B 110kg | 182cm

mmHg

*The pressure mapping scale is set to accurately perform with respect to the clinically recognised 32mmHg average capillary closing pressure. Clinicians should always ensure pressure mapping scale is set to accurately represent true performance.

Clinical Considerations

- Medical-grade foam for all day sitting
- No lateral movement in front castors for safe sit-to-stand transfer
- Tool-free adjustment
- High-level pressure-relieving cushion as standard
- Auto lockout of electronics when transport wheels are engaged



Resources
Clinical Assistance Guide
and Assessment Form

Performance Enhancing Accessories



Classic
Backrest



Advance
Immersion Backrest



Profiled
Headrest



Lateral
Support Blocks



Channel
Legrest

Petite	/	Black - CHP198293 Stone - CHP198294	Black - CHP198310 Stone - CHP198312	Black - CHP198315	Black - CHP198320
Small / Medium	Black - CHP198290	Black - CHP198295 Stone - CHP198297	Black - CHP198310 Stone - CHP198312	Black - CHP198315	Black - CHP198320
Medium / Large Tall	Black - CHP198291	Black - CHP198296 Stone - CHP198298	Black - CHP198311 Stone - CHP198313	Black - CHP198315	Black - CHP198320

Note: Blue highlights are used for identification only.

Technical Specifications

Code / Colour	Seat Height mm/”	Seat Depth mm/”	Seat Width mm/”	Backrest Height mm/”	Overall Height mm/”	Overall Depth mm/”	Overall Width mm/”	Distance from wall mm/”
Petite								
CHP198240 / Black CHP198241 / Stone	450	475 - 525	485 - 535	700	1080	900	760 - 810	580
CHP198240 - UK / Black CHP1982541 - UK / Stone	18	19 - 21	19 - 21	27.5	42.5	35.5	30 - 32	22.8
Small / Medium								
CHP198250 / Black CHP198251 / Stone	490	475 - 525	485 - 535	740	1140	900	760 - 810	580
CHP198250 - UK / Black CHP198251 - UK / Stone	19.3	19 - 21	19 - 21	29	45	35.5	30 - 32	22.8
Medium / Large Tall								
CHP198260 / Black CHP198261 / Stone	490	525 - 575	535 - 585	790	1180	950	820 - 870	580
CHP198260 - UK / Black CHP198261 - UK / Stone	19.3	21 - 23	21 - 23	31	46.5	37.4	32 - 34	22.8

Note: 60mm electric vertical height adjustment

BS 5852
PART 2 CRIB 5

2 YEAR WARRANTY
Frame & Electronics

1 YEAR WARRANTY
Handset & Cables



Australia
1300 133 120 | aspirecare.com.au

United Kingdom
01622 541235 | aidacare.co.uk