



aspire

# SEAT WALKER RANGE

aspire for...

Comfort  
Stability  
Support



[www.aspirecare.com.au](http://www.aspirecare.com.au)

1300 133 120



Aspire Seat Walkers provide assistance, support and relief; giving confidence to users of all shapes, sizes and mobility levels who want to remain actively mobile and independent.

### DELUXE SEAT WALKER



**Deluxe, In-Line Rear Wheels and Brakes**  
Increases space around the user's stepping zone

**Handle Height Range**  
835 to 990mm

**Taller Seat Height**  
590mm

**8" Solid Wheels**  
Ease over even the most robust indoor and outdoor terrain

Available colours: ■ WAF705750BL ■ WAF705750RE ✓ Deluxe is registered in the ARTG, AUST R ARTG 224650

### CLASSIC 8 SEAT WALKER



**Handle Height Range**  
820 to 975mm

**Taller Seat Height**  
575mm

**8" Solid Wheels**  
Ease over even the most robust indoor and outdoor terrain

Available colours: ■ WAF705700BL ■ WAF705700RE ✓ Classic 6 & 8 is registered in the ARTG, AUST R ARTG 224644

### CLASSIC 6 SEAT WALKER



**Handle Height Range**  
760 to 915mm

**Medium Seat Height**  
525mm

**6" Solid Wheels**  
Ease over every day surfaces indoors and out

Available colours: ■ WAF709600BL ■ WAF709600RE ✓ Classic 6 & 8 is registered in the ARTG, AUST R ARTG 224644

### MINI SEAT WALKER



**Handle Height Range**  
745 to 905mm

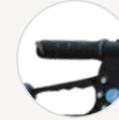
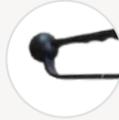
**Lower Seat Height**  
485mm

**6" Solid Wheels**  
Ease over every day surfaces indoors and out

Available colours: ■ WAF705800BL ■ WAF705800RE ✓ Mini is registered in the ARTG, AUST R ARTG 224651

## COMMON FEATURES & BENEFITS

✓ DVA Contracted ✓ AS / NZS Approved 1 YEAR LIMITED WARRANTY

									
<b>Rustproof Aluminium Frames</b> Are durable, stable and very lightweight	<b>Adjustable Handle Heights</b> Facilitate optimal sizing for postural positioning	<b>Latex-Free Ergonomically Shaped Grip Handles</b> Especially comfortable for users with arthritis	<b>Ultra-Soft Hand Brakes</b> With palm ball lock effortlessly	<b>Easy Folding Mechanism</b> Facilitates compact storage or transportation	<b>Padded Seats Flip Up</b> For access to clip on removable storage bags with discreet pocket compartments	<b>Curved Padded Backrests</b> For seated support and comfort	<b>360° Swivelling Front Wheels</b> For easy manoeuvrability	<b>Easy Wipe</b> Clean surfaces	<b>Tool Free Assembly and adjustment</b>

\*Common features listed on this spread are shared by the Classic 8", Classic 6", Mini, Deluxe, Push Down, XL and Superlite Adjustable Walkers

### PUSH DOWN SEAT WALKER



**Push Down Brakes**  
Are weight activated from sitting on the seat, ideal for arthritic hands

**Handle Height Range**  
855 to 965mm

**Seat Height**  
580mm

**7" Solid Wheels**  
Ease over every day surfaces indoors and out

Available colours: ■ WAF705600BL ■ WAF705600RE ✓ Push Down is registered in the ARTG, AUST R ARTG 224652

### SUPERLITE ADJUSTABLE SEAT WALKER



**Handle Height Range**  
740 to 920mm

**Adjustable Seat Height Range**  
490 to 595mm  
From low to tall

**6" Solid Wheels**  
Ease over every day surfaces indoors and out

**Lightest Seat Walker on the Market at only 5kg**

Available colour: ■ WAF709650CH ✓ Superlite Adj. is registered in the ARTG, AUST R ARTG 281703

### XL SEAT WALKER



**Handle Height Range**  
850 to 985mm

**Wide Padded Seat**  
With 180kg safe working load for large users

**Medium Seat Height**  
560mm

**8" Solid Wheels**  
Ease over even the most robust indoor and outdoor terrain

Available colour: ■ WAF750020SL ✓ XL is registered in the ARTG, AUST R ARTG 224645

### ACCESSORIES AVAILABLE

	
<b>OXYGEN BOTTLE HOLDER</b> Attaches to A-Frame Seat Walkers* WAA686905	<b>WALKING STICK HOLDER</b> Attaches To A-Frame Seat Walkers. Top mount clamp to securely hold the walking stick in place WAA686924
	
<b>SEAT WALKER BASKET</b> Clamp On Basket To Suit Aspire Seat Walkers WAA686910 - Deluxe WAA686909 - Classic, Mini & Push Down WAA686908 - Aspire XL	

\*Not compatible with XL Seat Walker (WAF750020SL)

## TRI WHEEL WALKER

The three wheel design makes the Tri Wheel Walker ideal for use in confined spaces. As well as a durable, robust steel construction, it also folds away for easy storage or transportation.



**Handle Height Range**  
825 to 925mm

**Convenient**  
Storage bag, basket and meal tray

**8" Solid Wheels**  
Provide smooth manoeuvrability over variable surfaces

Available colour:  
 WAF705500

 Tri Wheel is registered in the ARTG, AUST R ARTG 281703

## MEAL TRAY WALKER

Walking aid with meal trays, providing a means of transporting meals and other items within the home.



**Clip on Plastic Trays**  
With outer lip edge to capture spillage and prevent items from slipping off

**Height Adjustable Handles**  
With lockable handbrakes to hold the walker in place

**Large 6" Wheels**  
Make it easier to negotiate carpets and door thresholds

**360° Swivelling**  
Front wheels

Available colour:  
 DLG272061

 Meal Walker is registered in the ARTG, AUST R ARTG 271256

## MEAL TROLLEY

Meal transport trolley for use within the home, providing a means of transporting meals or other items from one room to another.



**Clip on Plastic Trays**  
With outer lip edge to capture spillage and prevent items from slipping off

**3" Castors**  
With two locking to hold the trolley in place

Available colour:  
 DLG272060

 Meal Trolley is registered in the ARTG, AUST R ARTG 271256

## COMMON FEATURES & BENEFITS:



**Easy Wipe**  
Clean surfaces



**Tool Free**  
Assembly and adjustment

 DVA Contracted

 AS / NZS Approved

**1 YEAR LIMITED WARRANTY**

## SPECIFICATIONS

### WALKERS

	OVERALL HEIGHT	OVERALL LENGTH	OVERALL WIDTH	SEAT DEPTH	SEAT WIDTH	SEAT HEIGHT	BACKREST HEIGHT	FOLDED WIDTH	WHEEL SIZE	UNIT WEIGHT	S.W.L
DELUXE	835 - 990mm	765mm	605mm	305mm	355mm	590mm	840mm	245mm	8"	7kg	130kg
CLASSIC 8	820 - 975mm	660mm	630mm	305mm	355mm	575mm	815mm	200mm	8"	7kg	130kg
CLASSIC 6	760-915mm	625mm	615mm	305mm	355mm	525mm	770mm	200mm	6"	6kg	130kg
MINI	745-905mm	560mm	615mm	305mm	355mm	485mm	690mm	200mm	6"	6kg	130kg
PUSH DOWN	855-965mm	660mm	620mm	305mm	355mm	580mm	815mm	255mm	7"	7kg	130kg
SUPERLITE ADJUSTABLE	740-920mm	635-765mm	500-670mm	305mm	355mm	490-595mm	700-870mm	280mm	6"	5kg	113kg
XL	850-985mm	700mm	735mm	340mm	495mm	560mm	850mm	230mm	8"	8kg	180kg
TRI WHEEL	825-925mm	575mm	715mm	N/A	N/A	N/A	N/A	275mm	8"	5.5kg	115kg

### MEAL WALKING AIDS

	OVERALL HEIGHT	OVERALL LENGTH	OVERALL WIDTH	LOWER TRAY HEIGHT	UPPER TRAY HEIGHT	TRAY WIDTH	WHEEL SIZE	UNIT WEIGHT	S.W.L
MEAL TRAY WALKER	770-920mm	560mm	585mm	415mm	770-920mm	450mm	6"	8kg	125kg
MEAL TROLLEY	810-940mm	505mm	460mm	355-480mm	615-740mm	450mm	3"	6kg	N/A

